
Instructions for use
EXPERTtorque E680 - 1.006.8700, 1.006.4700, 1.006.4300, 1.006.9000, 1.006.4600
EXPERTtorque Mini E677 - 1.007.3600, 1.006.0100, 1.006.3800

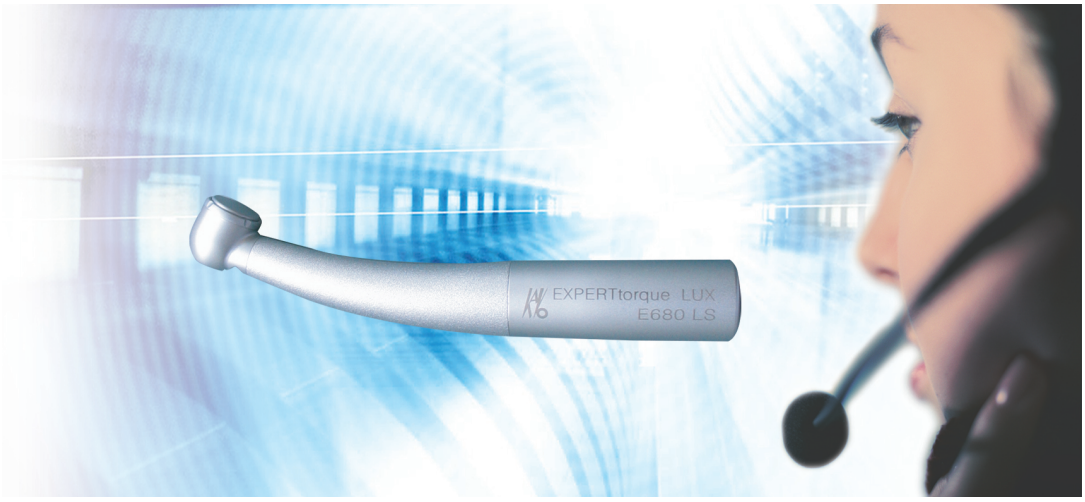

Always be on the safe side.

KaVo. Dental Excellence.

Distributed by:
KaVo Dental GmbH
Bismarckring 39
D-88400 Biberach
Phone +49 (0) 7351 56-0
Fax +49 (0) 7351 56-1488

Manufacturer:
Kaltenbach & Voigt GmbH
Bismarckring 39
D-88400 Biberach
www.kavo.com

CE 0124

Table of contents


1 User instructions 2
2 Safety 3
2.1 Description of safety instructions 3
2.2 Safety instructions 3
3 Product description 6
3.1 Purpose – Intended use 6
3.2 Technical Specifications 7
3.3 Transportation and storage conditions 7
4 Start up and shut down 9
4.1 Checking the amount of water 9
4.2 Checking the pressures 10
5 Operation 11
5.1 Attaching the medical device 11
5.2 Removing the medical device 11
5.3 Inserting the milling cutter or diamond grinder 11
5.4 Removing the milling tool or diamond grinder 12
6 Troubleshooting 13
6.1 Cleaning the spray nozzle 13
6.2 Changing the water filter 13
7 Reprocessing methods according to ISO 17664 14
7.1 Preparations at the site of use 14
7.2 Cleaning 14
7.2.1 Manual external cleaning 14
7.2.2 Automated external cleaning 14
7.2.3 Manual internal cleaning 15
7.2.4 Automated internal cleaning 15
7.3 Disinfection 15
7.3.1 Manual external disinfection 16
7.3.2 Manual internal disinfection 16
7.3.3 Automated external and internal disinfection 16
7.4 Drying 17
7.5 Care products and systems - Servicing 17
7.5.1 Servicing with KaVo Spray 17
7.5.2 Care with KaVo SPRAYrotor 18
7.5.3 Servicing with KaVo QUATTROcare 2104 / 2104A 18
7.5.4 Servicing with KaVo QUATTROcare PLUS 19
7.6 Packaging 19
7.7 Sterilisation 20
7.8 Storage 20
8 Tools and consumables 21
9 Terms and conditions of warranty 22

1 User instructions

Dear User
Congratulations on purchasing this KaVo quality product. By following the instructions below you will be able to work smoothly, economically and safely.

© Copyright by KaVo Dental GmbH

Symbols

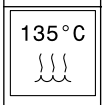	Refer to the chapter on Safety/Warning symbol
	Important information for users and service technicians
	Can be steam-sterilised at 134 °C -1 °C / +4 °C (273 °F -1.6 °F / +7.4 °F)
	Thermodisinfectable
	CE mark (European Community). A product bearing this mark meets the requirements of the applicable EC directive.
▶	Action request

Target group

This document is intended for dentists and their assistants. The section on starting up is also intended for service technicians.

2 Safety

2.1 Description of safety instructions

	Warning symbol

Structure

	⚠ DANGER **The introduction describes the type and source of the hazard.** This section describes potential consequences of non-compliance. ▸ The optional step includes necessary measures for hazard prevention.

Description of hazard levels

The safety instructions listed here, together with the three levels of danger will help avert property damage and injury.

	⚠ CAUTION CAUTION indicates a hazardous situation that can cause damage to property or mild to moderate injuries.

	⚠ WARNING WARNING indicates a hazardous situation that can lead to serious or fatal injury.

	⚠ DANGER DANGER indicates a maximal hazard due to a situation that can directly cause death or fatal injury.

2.2 Safety instructions

	⚠ WARNING **Hazards for the care provider and the patient.** In the case of damage, irregular running noise, excessive vibration, un-typical warming or when the cutter or grinder cannot be held. ▸ Do not use further and notify Service.



⚠ WARNING

Hazard from contraindication.
If the soft tissue in the oral cavity is injured, the compressed air may enable septic substances to enter the tissue.

- ▸ Treatment must be discontinued with instruments operated by compressed air when soft tissue is damaged in the oral cavity.



⚠ CAUTION

Risk due to incorrectly stored instrument.
Injury and infection caused by chucked cutters or grinders.
Damage to clamping system from dropping the instrument.

- ▸ After treatment, place the instrument properly in the cradle, without the cutter or grinder.



⚠ CAUTION

Risk of burn injury from hot handpiece head or hot handpiece lid.
Burn injuries in the mouth may be caused if the handpiece overheats.

- ▸ Never touch soft tissue with the handpiece head or handpiece lid!



⚠ CAUTION

Hazard from use as a light probe.
Do not use the device as a light probe since the rotating cutters or grinders can cause injury.

- ▸ Use an appropriate light probe for additional illumination of the oral cavity or site of preparation.



⚠ CAUTION

Premature wear and malfunctioning from improper storage during long periods of nonuse.
Reduced product life.

- ▸ The medical device should be cleaned, serviced and stored in a dry location, according to instructions, before long periods of nonuse.



Note


For safety reasons, we recommend that the tool holder system be checked annually after the warranty period expires.

The following individuals are authorized to repair and service KaVo products:
- Technicians at KaVo branches throughout the world
- Technicians specially trained by KaVo

To ensure proper function, the medical device must be set up according to the reprocessing methods described in the KaVo Instructions for Use, and the care products and care systems described therein must be used. KaVo recommends specifying a service interval at the dental office for a licensed shop to clean, serv-


ice and check the functioning of the medical device. This service interval depends on the frequency of use and should be adjusted accordingly.
Service may only be carried out by KaVo-trained repair shops using original KaVo replacement parts.

3 Product description

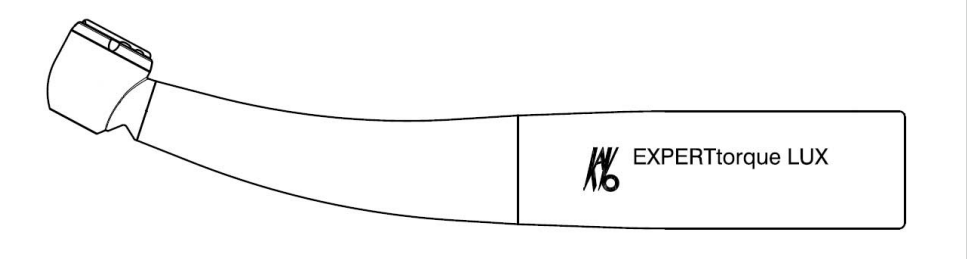

- EXPERTtorque LUX E680 L**Mat. no. 1.006.8700**
- EXPERTtorque E680 C**Mat. no. 1.006.9000**

- EXPERTtorque LUX E680 LN **Mat. no. 1.006.4700**
- EXPERTtorque LUX E680 LM **Mat. no. 1.006.4300**
- EXPERTtorque LUX E680 LS **Mat. no. 1.006.4600**

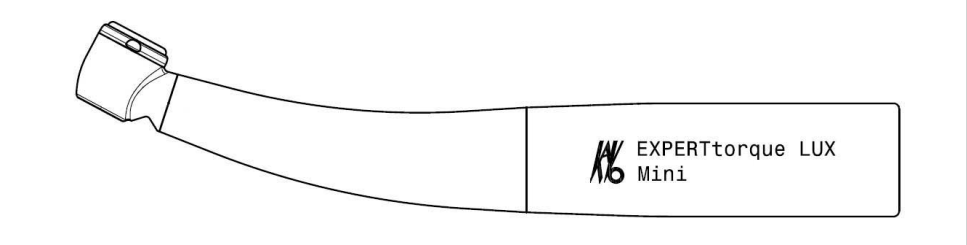

- EXPERTtorque Mini LUX E677 L**Mat. no. 1.007.3600**

- EXPERTtorque Mini LUX E677 LN **Mat. no. 1.006.0100**
- EXPERTtorque Mini LUX E677 LM **Mat. no. 1.006.3800**

3.1 Purpose – Intended use

Purpose:

This medical device is
- Only intended for dental treatment. Any other type of use or alteration to the product is impermissible and can be hazardous. The medical device is intended for the following uses: Removal of carious material, cavities and crown preparations, removal of fillings, processing of tooth, restoration surfaces and surgical procedure.
- A medical device according to relevant national statutory regulations.

Proper use:

According to these regulations, this medical device may only be used for the described application by a knowledgeable user. The following must be observed:
- the applicable health and safety regulations
- the applicable accident prevention regulations
- these instructions for use

According to these regulations, it is the responsibility of the user to:
- only use equipment that is operating correctly,
- use the equipment for the proper purpose,
- protect him or herself, the patient and third parties from danger, and
- avoid contamination from the product.

3.2 Technical Specifications

	KaVo MULTIflex	NSK PTL-CL-LED3	Morita Alpha Connection	Sirona R/F
Drive pressure (bar)	2.1 - 3.5	2.0 - 3.0	3.5 - 3.7	2.6 - 3.0
Drive pressure recommended (bar)	> 2.8	2.5	3.6	2.7
Return air pressure (bar)	< 0.5	< 0.5	< 0.3	< 0.3
Spray water pressure (bar)	0.8 - 2.0	0.8 - 2.0	0.5 - 2.0	0.8 - 2.2
Spray air pressure (bar)	1.0 - 2.5	1.0 - 2.5	-	0.8 – 2.9

	KaVo MULTIflex	NSK PTL-CL-LED3	Morita Alpha Connection	Sirona R/F
Air consumption (Nl/min)	39 - 51			
Idling speed (rpm) EXPERTtorque / EXPERTtorque Mini	340,000 - 420,000 / 400,000- 480,000			
Recommended application force (N)	2 - 3			

Note

KaVo MULTIflex is a registered trademark of Kaltenbach & Voigt GmbH, Biberach.
NSK PTL - CL- LED 3 is a product name of NSK NAKANISHI INC., Japan.
Morita Alpha Connection is a product name of MORITA Dental Company, Japan.
Sirona is a registered brand of Sirona Dental System GmbH, Bensheim.

3.3 Transportation and storage conditions

⚠ CAUTION

It is hazardous to start up the medical device after it has been stored strongly refrigerated.
This can cause the medical device to malfunction.

- ▸ Prior to start-up, very cold products must be heated to a temperature of 20°C to 25°C (68°F to 77°F).

	Temperature: -20°C to +70°C (-4°F to +158°F)
	Relative humidity: 5% RH to 95% RH absence of condensation
	Air pressure: 700 hPa to 1060 hPa (10 psi to 15 psi)
	Protect from moisture

4 Start up and shut down



⚠ WARNING

Hazard from non-sterile products.
Infection hazard for care provider and patient.

- Before first use and after each use, prepare and sterilise the medical device and accessories accordingly.



⚠ WARNING

Disposal of the product in the appropriate manner.
Prior to disposal, the product and accessories must be appropriately prepared or sterilised if this is necessary.



⚠ CAUTION

Damage from contaminated and moist cooling air.
Contaminated and moist cooling air can cause malfunctions.

- Make sure that the supply of cooling air is dry, clean, and uncontaminated according to EN ISO 7494-2.


4.1 Checking the amount of water



⚠ CAUTION

Overheating of the tooth due to lack of cooling water.
Insufficient spray water can cause the medical device to overheat and damage the pulp and tooth.

- Adjust the water amount for the spray head to a min. of 50 cm^3/min!
- Check the spray water channels and clean the spray nozzles with the nozzle needle **Mat. no. 0.410.0921** if necessary.

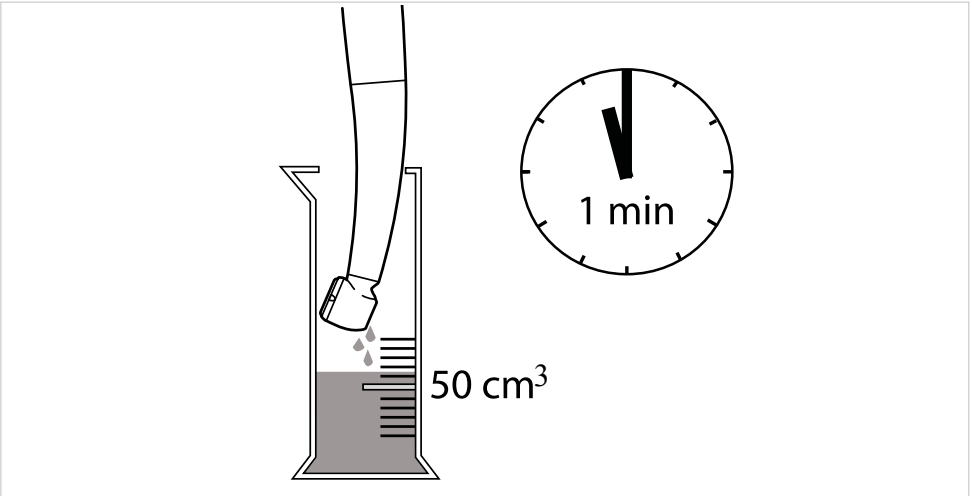

4.2 Checking the pressures



⚠ CAUTION

Compressed air connection on devices.
Contaminated and humid compressed air leads to premature wear and tear.

- Ensure that the cooling air is dry, clean and uncontaminated in accordance with EN ISO 7494-2.



- Insert the test manometer (**Mat. no. 0.411.8731**) between the coupling and the medical device and check the following pressures:

  Drive pressure, drive pressure recommended, return air pressure, spray water pressure, and spray air pressure.

See also: 3.2 Technical Specifications, Page 7

A higher drive pressure will be reduced automatically by the medical device.

5 Operation

Note

At the beginning of each workday, the water-conducting systems should be rinsed for at least 2 min. without the instrument being attached; if there is a risk of contamination from reflux or back suction, the system must be rinsed for 20-30 seconds.

5.1 Attaching the medical device

⚠ WARNING

Detachment of the medical device during treatment.
A medical device that is not properly locked can detach from the coupling during treatment.

- Before each use, check if the medical device is securely locked onto the coupling.

⚠ CAUTION

Damage from inaccurate coupling.
Inaccurate coupling (especially during the afterglow period) can destroy the high-pressure lamp or the LED of the coupling or reduce its service life.

- Check the seat of the turbine on the coupling by pulling on it.



- Accurately attach the medical device to the coupling and push is to the rear until the coupling audibly locks in the medical device.
- Pull on it to make sure that the medical device is securely affixed to the coupling.

5.2 Removing the medical device

- Hold the coupling tight, and pull the medical device off while twisting slightly.

5.3 Inserting the milling cutter or diamond grinder

Note

Only use carbide cutters or diamond grinders that comply with EN ISO 1797-1 type 3, are made of steel or hard metal and meet the following criteria:
- Shaft diameter: 1.59 to 1.60 mm
- Overall length: EXPERTtorque: max. 25 mm
- Overall length of EXPERTtorque Mini: max. 19 mm
- Shaft clamping length: EXPERTtorque: min. 11 mm
- Shaft clamping length: EXPERTtorque Mini: min. 9 mm
- Blade diameter: max. 2 mm



⚠ WARNING

Use of unauthorised cutters or grinders.
Injury to the patient or damage to the medical device.

- Observe the instructions for use and use the cutter or grinder properly.
- Only use cutters or grinders that do not deviate from the specified data.



⚠ CAUTION

Injury from using worn drill bits or burs.
Drill bits or burs could fall out during treatment and injure the patient.

- Never use drill bits or burs with worn shafts.



⚠ CAUTION

Danger of injury from cutters or grinders.
Infections or cuts.


- Wear gloves or fingerstalls.



⚠ CAUTION

Hazard from defective chucking system.
The cutter or grinder could fall out and cause injury.

- Pull on the cutter or grinder to check that the chucking system is okay and the cutter or grinder is securely held. When checking, inserting and removing, use gloves or a fingerstall to prevent an injury or infection.

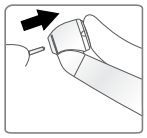

- Forcefully press the push button with your thumb and simultaneously insert the cutter or grinder all the way.
- Check that the cutter or grinder is seated by pulling on it.


5.4 Removing the milling tool or diamond grinder



⚠ WARNING

Hazard from rotating cutter or grinder.
Lacerations and damage to the chucking system.

- Do not touch rotating cutter or grinder.!
- Never press the press-button while the cutter or grinder is rotating!
- Remove the cutter or grinder from the contra-angle handpiece after treatment to avoid injury or infection while storing it.

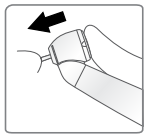

- After the cutter or grinder has stopped rotating, press the press-button with your thumb and simultaneously pull out the drill bit or bur.

6 Troubleshooting


6.1 Cleaning the spray nozzle



⚠ CAUTION


Hazard from insufficient spray water.
Overheating of the medical device and damage to the tooth.

- Check the spray water channels and clean the spray nozzles with the nozzle needle **Mat. no. 0.410.0921** if necessary.

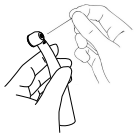


- Use the nozzle needle (**Mat. no. 0.410.0921**) to free the water passage in the spray nozzles.

6.2 Changing the water filter

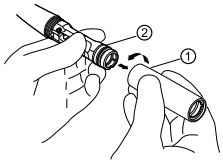

- Unscrew the sleeve ① in counterclockwise direction from the insert ② and pull it off.

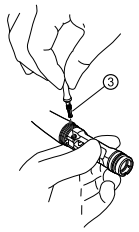

- Unscrew the water filter ③ with the wrench (Mat. no. 1.002.0321) and take it out

- Insert the new filter (Mat. no. 1.002.0271) and screw it in with the wrench
- Place the sleeve ① on the insert ②, and screw it tight in clockwise direction.

7 Reprocessing methods according to ISO 17664

7.1 Preparations at the site of use



⚠ WARNING

Hazard from nonsterile products.
There is a risk of infection from contaminated medical devices.

- ▸ Take suitable personal protective measures.

- ▸ Remove all residual cement, composite or blood without delay.
- ▸ Recondition the medical device as soon as possible after treatment.
- ▸ Remove the cutter or grinder from the medical device.
- ▸ The medical device must be dry when transported for reconditioning.
- ▸ Do not place it in a solution or similar.

7.2 Cleaning




⚠ CAUTION

Malfunctions from cleaning in the ultrasonic unit.
Defects in the product.

- ▸ Only clean manually or in a thermodisinfector.

7.2.1 Manual external cleaning

Accessories required:
- Tap water 30 °C ± 5 °C (86 °F ± 10 °F)
- Brush, e.g. medium-hard toothbrush

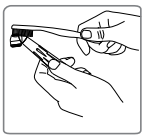

- ▸ Brush off under flowing tap water.

7.2.2 Automated external cleaning



KaVo recommends thermodisinfectors in accordance with EN ISO 15883-1, which are operated with alkaline cleaning agents with a pH value of max. 10 (e.g. Miele G 7781/G 7881 – Validation was carried out with Programme "VARIO-TD", cleaning agent "neodisher® mediclean", neutralisation agent "neodisher® Z" and rinsing agent "neodisher® mielclear" and only applies to the material compatibility with KaVo products).

- ▸ For program settings as well as cleansers and disinfectants to be used, please refer to the Instructions for Use of the thermodisinfector.

- ▸ Directly after automated cleaning, treat the medical device with the care products and systems provided by KaVo.

7.2.3 Manual internal cleaning

Can only be done with KaVo CLEANspray or KaVo DRYspray.

Currently, this is feasible only with KaVo turbines equipped with the MULTIflex coupling.

- ▸ Cover the medical device with the KaVo CLEANpac bag, and place it on the corresponding care adapter. Press the spray button three times for 2 seconds each time. Remove the medical device from the spray attachment and let the cleaner work for one minute.
- ▸ Afterwards, rinse for 3-5 seconds with KaVo DRYspray.

See also: KaVo CLEANspray / KaVo DRYspray Instructions for Use

Note

KaVo CLEANspray and KaVo DRYspray for manual interior cleaning are only available in the following countries:
Germany, Austria, Switzerland, Italy, Spain, Portugal, France, Luxembourg, Belgium, Netherlands, United Kingdom, Denmark, Sweden, Finland and Norway.
In other countries interior cleaning can only be carried out with thermodisinfectors in accordance with EN ISO 15883-1.

7.2.4 Automated internal cleaning

KaVo recommends thermodisinfectors in accordance with EN ISO 15883-1, which are operated with alkaline cleaning agents with a pH value of max. 10 (e.g. Miele G 7781/G 7881 – Validation was carried out with Programme "VARIO-TD", cleaning agent "neodisher® mediclean", neutralisation agent "neodisher® Z" and rinsing agent "neodisher® mielclear" and only applies to the material compatibility with KaVo products).

- ▸ For program settings as well as cleansers and disinfectants to be used, please refer to the Instructions for Use of the thermodisinfector (complying with max. pH value of 10).

- ▸ In order to prevent negative effects on the medical device, make sure that the interior and the exterior of the medical device are dry after completion of the cycle, and then grease it immediately with servicing agents from the KaVo care system.


7.3 Disinfection

⚠ CAUTION

Malfunctioning from using a disinfectant bath or disinfectant containing chlorine.
Defects in the product.

- ▸ Only disinfect in a thermodisinfector or manually.

7.3.1 Manual external disinfection

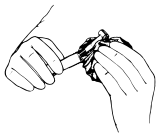

KaVo recommends the following products based on material compatibility.

- Mikrozid AF made by Schülke & Mayr (liquid or cloths)
- FD 322 made by Dürr

- CaviCide made by Metrex

Consumables required:
- Cloths for wiping off the medical device.

- ▸ Spray the disinfectant on a cloth, then thoroughly wipe down the medical device and leave the disinfectant to soak in according to the instructions from the disinfectant manufacturer.

- ▸ Follow the instructions for use of the disinfectant.

7.3.2 Manual internal disinfection

The efficacy of manual internal disinfection must be demonstrated by the manufacturer of the disinfection agent. With KaVo products, use only disinfection agents that have been released by KaVo with respect to the compatibility of materials (e.g. WL-cid / made by ALPRO).

- ▸ Follow the instructions for use of the disinfectant.

- ▸ Immediately after internal disinfection, lubricate the KaVo medical device immediately with care agents from the KaVo care system.

7.3.3 Automated external and internal disinfection



KaVo recommends thermodisinfectors in accordance with EN ISO 15883-1, which are operated with alkaline cleaning agents with a pH value of max. 10 (e.g. Miele G 7781/G 7881 – Validation was carried out with Programme "VARIO-TD", cleaning agent "neodisher® mediclean", neutralisation agent "neodisher® Z" and rinsing agent "neodisher® mielclear" and only applies to the material compatibility with KaVo products).

- ▸ For program settings as well as cleansers and disinfectants to be used, please refer to the Instructions for Use of the thermodisinfector (complying with max. pH value of 10).

- ▸ In order to prevent negative effects on the medical device, make sure that the interior and the exterior of the medical device are dry after completion of the cycle, and then grease it immediately with servicing agents from the KaVo care system.

7.4 Drying

Manual Drying

- Blow off the outside and inside with compressed air until water drops are no longer visible.

Automatic Drying

The drying procedure is normally part of the cleaning program of the thermodisinfector.

Note
Please observe the instructions for use of the thermodisinfector.

7.5 Care products and systems - Servicing


⚠ WARNING
Sharp cutters or grinders in the medical device.
Risk of injury from sharp or pointed cutters or grinders.
- Remove cutter or grinder.

⚠ CAUTION
Premature wear and malfunctions from improper servicing and care.
Reduced product life.
- Perform regular proper care and maintenance!

Note
KaVo only guarantees that its products will function properly when the care products used are those listed as accessories, as they were tested for proper use on our products.

7.5.1 Servicing with KaVo Spray


KaVo recommends servicing the product after each time it is used, i.e. after each automatic cleaning and before each sterilisation.

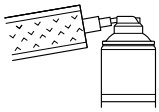

- Remove the cutter or grinder.
- Cover the product with the Cleanpac bag.
- Plug the product onto the cannula, and press the spray button for one second.

Chuck care


KaVo recommends cleaning and servicing the chuck system once a week.

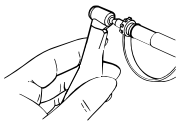

- Remove the cutter or grinder, place the spray nipple tip in the opening and spray.

- Carry out the servicing according to the instructions in the section, "Care with KaVo Spray".

7.5.2 Care with KaVo SPRAYrotor

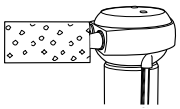


KaVo recommends servicing the product after each time it is used, i.e. after each automatic cleaning and before each sterilisation.

- Place the product on the appropriate coupling of the **KaVo SPRAYrotor** and cover it with a CLEANpac bag.
- Service the product.

See also: Instructions for use **KaVo SPRAYrotor**

7.5.3 Servicing with KaVo QUATTROcare 2104 / 2104A

Cleaning and care unit with expansion pressure for effective cleaning and care.

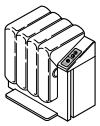

KaVo recommends servicing the product as part of the reprocessing after each use, i.e. after each cleaning, disinfection, and before each sterilisation.

- Remove the cutter or grinder.
- Servicing the product.

Chuck care

KaVo recommends cleaning and servicing the chuck system once a week.

See also: Instructions for use **KaVo QUATTROcare 2104 / 2104A**




- Remove the cutter or grinder, place the spray nipple tip in the opening and spray.

- Subsequently treat with the care products and care systems specified.

See also: Servicing with KaVo **QUATTROcare 2104 / 2104A**

7.5.4 Servicing with KaVo QUATTROcare PLUS

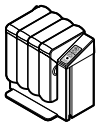

KaVo recommends servicing the product after each time it is used, i.e. after each automatic cleaning and before each sterilisation.

- Remove the cutter or grinder.
- Servicing the product in QUATTROcare PLUS.

See also: Instructions for Use **KaVo QUATTROcare PLUS 2124 A**

Servicing of the clamping chuck


KaVo recommends cleaning and servicing the chuck system once a week.

See also: Instructions for Use **KaVo QUATTROcare PLUS 2124 A**

Note

Handpieces must be taken off the service couplings before the chuck service can be started and performed.

- Remove the service coupling chuck from the side hatch of the QUATTROcare PLUS and attach it to coupling service point four, on the far right. A MULTIflex adapter must be mounted there.

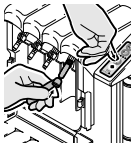

- Press the handpiece together with the guide bush of the chuck to be serviced against the tip of the service coupling chuck.
- Press the button showing the chuck service symbol.

Note

Close the chuck service mode.
Option 1: Place the dental handpieces in the QUATTROcare PLUS 2124 A, close the front door and start theservice procedure.
Option 2: After three minutes with no service procedure running, the device automatically switches back to normal service mode.

See also: Care with KaVo QUATTROcare PLUS

7.6 Packaging

Note

The sterilisation bag must be large enough for the handpiece so that the bag is not stretched.
The quality and use of the sterilisation packaging must satisfy applicable standards and be suitable for the sterilisation procedure!

- Individually seal the medical device in the sterilised item packaging.

7.7 Sterilisation

Sterilisation in a steam steriliser (autoclave) in accordance with EN 13060 / ISO 17665-1



⚠ CAUTION

Premature wear and malfunctions from improper servicing and care.
Reduced product life.

- ▸ Before each sterilisation cycle, service the medical device with KaVo care products.



⚠ CAUTION

Contact corrosion due to moisture.
Damage to product.

- ▸ Immediately remove the product from the steam steriliser after the sterilisation cycle!



The medical device has a maximum temperature resistance up to 138 °C (280.4 °F).

Select a suitable procedure (depending on the available autoclave) from the following sterilisation processes:

- Autoclave with three times pre-vacuum:
  - at least 3 minutes at 134 °C -1 °C / +4 °C (273 °F -1.6 °F / +7.4 °F)
- Autoclave using the gravity method:
  - at least 10 minutes at 134 °C -1 °C / +4 °C (273 °F -1.6 °F / +7.4 °F)

- ▸ Use according to the manufacturer's Instructions for Use.

7.8 Storage

- ▸ Reprocessed products should be stored protected from dust with minimum exposure to germs in a dry, dark and cool space.

- ▸ Comply with the expiry date of the sterilised items.

8 Tools and consumables

Available from dental suppliers.

Material summary	Mat. no.
Replacement turbine EXPERTtorque	1.007.9313
Replacement turbine EXPERTtorque Mini	1.007.9457
Wrench for lid of EXPERTtorque	0.411.3053
Wrench for lid of EXPERTtorque Mini (new)	1.008.6133
Wrench for lid of EXPERTtorque Mini (old)	1.006.3384
Replacement filter	1.002.0271
Wrench	1.002.0321
Instrument stand 2151	0.411.9501
Insert for turbines	0.411.9902
Nozzle pin	0.410.0921
Cleanpac 10 units	0.411.9691
KaVo MULTIflex spray head for KaVo Spray	0.411.9921
NSK spray head for KaVo Spray	1.005.8436
Morita spray head for KaVo Spray	1.005.8275
SIRONA spray head for KaVo Spray	1.005.8365
KaVo MULTIflex service coupling for KaVo QUATTROcare	0.411.7991
NSK service coupling for KaVo QUATTROcare	1.000.8786
Morita service coupling for KaVo QUATTROcare	1.000.6063
Sirona service coupling for KaVo QUATTROcare	1.000.7156

Material summary	Mat. no.
Adapter KaVo MULTIflex for KaVo CLEANspray/ DRYspray	1.007.1775
CLEANspray/ DRYspray Starter set 2116 P	1.007.0573
KaVo CLEANspray 2110 P	1.007.0579
KaVo DRYspray 2117 P	1.007.0580
KaVo Spray 2112 A	0.411.9640
ROTAspray 2142 A	0.411.7520
QUATTROcare plus Spray 2140 P	1.005.4525

9 Terms and conditions of warranty

The following warranty conditions apply to this KaVo medical device:

KaVo provides the end customer with a warranty of proper function and guarantees zero defects in respect of material and processing for a period of 18 months from the date of the invoice, subject to the following conditions:
In case of justified complaints, KaVo will honour its warranty with a free replacement or repair. Other claims of any nature whatsoever, in particular with respect to compensation, are excluded. In the event of default, gross negligence or intent, this shall only apply in the absence of mandatory legal regulations to the contrary.
KaVo shall not be liable for defects and their consequences that have arisen or may arise from natural wear, improper handling, cleaning or maintenance, non-compliance with operating, maintenance or connection instructions, calcination or corrosion, contaminated air or water supplies or chemical or electrical factors deemed abnormal or impermissible in accordance with KaVo's instructions for use or other manufacturer's instructions. The warranty granted does not usually extend to lamps, light conductors made of glass and glass fibres, glassware, rubber parts, and the colourfastness of plastic parts.
All liability is excluded if defects or their consequences originate from manipulations or changes to the product made by the customer or a third party that is not authorised by KaVo.
Warranty claims will only be accepted if the product is submitted along with proof of purchase in the form of a copy of the invoice or note of delivery. The dealer, purchase date, type, and serial number must be clearly evident from this document.

1.007.9662 · kb · 20130123 - 05 · en

KaVo. Dental Excellence.