ULTRADENT PRODUCTS, INC.

PQ1®

Single Component Bonding Agent

Description:
The PQ1® single syringe delivered, single component resin bonding system.
1. PQ1 is a single syringe delivered, light cured bonding agent with ethyl alcohol solvent carrier. It will cure with most all lights (not laser). PQ1 is 40% filled and is radiopaque.
Note: Upon discharge from the syringe PQ1 will appear milky; after proper air thinning, PQ1 becomes translucent.
2. Ultra-Etch® 35% phosphoric acid is unsurpassed for dentin and enamel conditioning and is self-limiting in its depth of etch (averge of 1.9µm).

Indications:
Use for most bonding needs in restorative dentistry. PQ1 is formulated to work well with flowable composite for light cure luting.

PQ1 bonds to the following materials:
1. Dentin and enamel, section **I**.
2. Porcelain, section **II**.
3. Metal, section **III**.
4. Composite, section **IV**.

Fig. 1 Etch: *Apply Ultra-Etch 20 seconds to dentin/enamel.Wash with air/water spray. Blot or blow off all surface water—do not desiccate.*

Fig. 2 Bond: *Apply PQ1 primer/ bonding resin to dentin and enamel for 10 seconds. Air thin 10 seconds to the point of losing milky appearance,* ***not too thin.***

Fig. 3 Restore: *Light cure 20 seconds if using a standard halogen light, or 10 seconds if using a high energy light. Restore with quality composite.*

Procedure:
Remove luer lock caps from syringes and attach appropriate tips: ***Ultra-Etch:*** Inspiral® Brush tip (some prefer the Blue Micro® tip). ***PQ1:*** Inspiral Brush tip or Black Micro® FX™ tip.

I. Dentin and Enamel
1. ***Prepare Surface:*** Prepared dentin and enamel surfaces to be bonded should be clean and isolated from contamination from soft tissue and oral fluids. Make sure surfaces are caries free. For abrasion/abfraction Class V preparations, roughen with a diamond bur.
2. ***Etch:*** Apply Ultra-Etch onto preparation surfaces and leave for 20 seconds (Fig. 1). Rinse dentin and enamel thoroughly for 5 seconds with firm air/water spray. With brief air blasts, remove all visible water. Do not desiccate.
3. ***Bond:*** Immediately apply a uniform coating of PQ1 to the surface (Fig. 2). Rub with moderate pressure on the dentin for 10 seconds. Avoid rubbing on enamel. Thin/dry with a focused air stream with 1/2 air pressure at 10mm for 10 seconds to remove alcohol solvent. Surface should look glossy but not puddled or milky.
Note: An alternative to air thinning is to vacuum thin with Ultradent's luer vacuum adapter and an 18-22 gauge Endo-Eze® cannula. This prevents over thinning and splatter of resin to soft tissues. Follow with 5-10 seconds of gentle air drying to evaporate solvent carrier.

4. *Light Cure:* Light cure 20 seconds if using a quality standard halogen light. Light cure 10 seconds if using a high energy light such as VALO. *Note: PO1 has been formulated to cure with most light sources. Air thin soon after applying. Prolonged exposure to ambient light may cause resin to polymerize prematurely. The operatory light may need to be turned down. We recommend covering syringe tip if left exposed to ambient light for long periods of time during a procedure.*

A. *For Composite:* Flowable composite may be used for a superadaptive initial layer. Follow with paste composite. Use a paste composite for the last 2-3mm layer. For surfaces subject to minimal wear such as Class V's, a flowable composite may be used to fill the surface.

B. *For Amalgam:* Follow steps 1, 2, and 3 then simply condense mixed amalgam into cavity preparation. Oxygen inhibition layer provides mechanism for bond since polymerization occurs under anaerobic conditions as amalgam fills the preparation. For highest bond strength with fresh placed amalgam, we recommend PQ Amalgam.

II. Porcelain Repair (Figs. 4-6):
Remove luer lock caps and attach the appropriate tips: ***Ultradent® Porcelain Etch:*** Inspiral Brush tip. ***Ultradent® Silane:*** Black Mini Brush tip. Before clinical procedure, dispense Ultradent Porcelain Etch onto paper pad to verify evenness of flow. If there is resistance to flow, DO NOT proceed. Replace with new tip and verify again before clinical use.
1. *Prepare Surface:* Remove loose and/or thin projections of porcelain and bevel as needed. Freshen fractured surface of porcelain with medium diamond bur and/or microabrade.
Note: Ultradent Porcelain Etch technique has been shown to be superior to microabrasion with feldspathic porcelains.
2. *Isolate:* Protect adjacent tissues and restorations before etching. Use OraSeal® Caulking to ensure rubber dam seal. When rubber dam placement is not possible, isolation can be established with EtchArrest acid barrier/neutralizer (Fig. 5). Also use EtchArrest to neutralize hydrofluoric acid before rinsing to ensure safe, controlled removal.
3. *Etch:* Apply Porcelain Etch with the Inspiral Brush tip (Figs. 4 and 5) to prepared surface. Etch 90 seconds, then neutralize with EtchArrest. Remove material with vacuum and gentle rinse, then 5 second air/water wash, then dry.
4. *Silanate:* Brush Silane into etched porcelain. Let stand for 60 seconds and air dry to evaporate any remaining carrier.
5. *Bonding Agent:* Apply a puddled coat of PQ1 to silanated surface. Thin/dry with a focused air stream with 1/2 air pressure at 10mm for 10 seconds to remove alcohol solvent. Surface should look glossy but not puddled or milky. Light cure 20 seconds if using a standard halogen light, 10 seconds if using a high energy light.
6. *Restore:* Use appropriate composite, building incrementally.

III. Metal to cover/repair:
1. *Prepare Surface:* Microabrade metal surface. Note Fig. 6. *Note: Bond strength to metal is determined by metal alloy and surface preparation technique.*
2. *Clean:* Metal surface should be cleaned with Ultra-Etch for 10 seconds followed by a 5 second rinse with firm air/water spray. Dry thoroughly.
3. *Bonding Agent:* Apply a puddled coat of PQ1 to etched surface for 10 seconds. Thin/dry with a focused air stream with 1/2 air pressure at 10mm for 10 seconds to remove alcohol solvent. Surface should look glossy but not puddled or milky. Light cure 20 seconds if using a standard curing light, or 10 if using a high energy curing light such as VALO.
4. *Restore:* Composite resin of choice may now be placed. Start with 0.5-1.5mm initial layer. Cure, then build in incremental layers 2-2mm thick.
Note: When metal opaquing or opaquing of dark dentin is necessary, use a thin layer of PermaFlo® Dentin Opaquer and light cure before placement of composite resin.

IV. Composite Repair:
1. *Prepare Surface:* Freshen with diamond bur. With older macro-filled composites, porcelain etching for 90 seconds or microabrasion increases bond strength. Remove weakened portions of existing composite. Open defect to limit risk of air entrapment during composite placement (using microabrasion provides superior results).
2. *Etch/Clean:* The prepped composite surface should be cleaned with Ultra-Etch for 10 seconds followed by a 5 second rinse with firm air/water spray. Dry thoroughly.
3. *Bonding Agent:* Apply a puddled coat of PQ1 to entire prepped area. Rub with moderate pressure for 10 seconds. Thin/dry with a focused air stream with 1/2 air pressure at 10mm for 10

Fig. 4 *Isolate with a rubber dam whenever possible.*

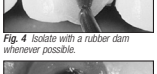
Fig. 5 *EtchArrest should be placed when rubber dam placement is not practical. EtchArrest is placed to protect soft tissues and adjacent porcelain. EtchArrest neutralizes hydrofluoric acid upon contact.*

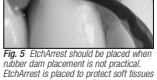
Fig. 6 *Here porcelain is etched for maximum bond strength to porcelain.*

seconds to remove alcohol solvent. Surface should look glossy, not puddle or milky. Light cure 20 seconds if using a standard halogen light, 10 seconds if using a high energy light such as Valo.
4. ***Restore:*** Build incrementally with an appropriate composite.

V. Light Cure Luting:
1. ***Prepare surface:*** Scour preparations with Consepsis® Scrub or pumice using Star brush. Retract tissue and control flow of sulcular fluid. Make sure porcelain veneer is clean, etched and silanated.
2. ***Etch:*** 20 Second Etch on enamel and dentin with Ultra-Etch. Rinse firmly with air/water spray 5 seconds. Blot or blow all visible water off but do not desiccate.
3. ***Bonding Agent:*** Apply a puddled coat of PQ1. Thin/dry with a focused air stream with 1/2 air pressure at 10mm for 10 seconds to remove alcohol solvent. Surface should look glossy but not puddled or milky. DO NOT light cure.
4. ***Lute:*** Apply luting/restorative flowable composite to silanated veneer; seat to place, wipe off excess resin. Apply DeOx to prevent oxygen inhibition on margins then light cure. PQ1 is polymerized simultaneously with luting/ composite resin.

Precautions:
A. Ultra-Etch & Porcelain Etch:
1. Remember to establish adequate isolation (DermaDam, cotton rolls, etc.) before using etchants. Isolation is also important to prevent contamination during bonding agent application. Use EtchArrest to protect tissues and neutralize hydrofluoric acid for situations where rubber dam is not feasible. Any tissues or surfaces inadvertently exposed to porcelain etch should be immediately rinsed with copious water.
2. It is recommended that the clinician verify syringe flow from tip on a mixing pad BEFORE using on patients. Before application, always verify tip is firmly attached to the luer lock assembly.
3. The Inspiral Brush or Blue Micro tip are designed and recommended for application of Ultradent etchants. Use of inappropriate tips may cause extreme flow resistance. If ease of flow is inhibited, replace tip. Forcing syringe plunger could cause uncontrolled discharge of syringe contents.
4. To achieve quality dentin bonding, do not desiccate dentin prior to PQ1 application. Remove excess water with controlled air blasts, by blotting, or with the Luer Vacuum adapter using a 18-22 gauge Endo-Eze cannula.
B. PQ1 Bonding Agent:
1. Resins can be sensitizing. Dental personnel should avoid repeated contact. Do not use on patients with a known sensitivity to acrylates or other resins.
2. Light-activated bonding agents are sensitive to ambient light; syringe caps should be replaced following use to prevent unwanted polymerization. We recommend covering syringe tip with 2 x 2 gauze if left exposed to ambient light for long periods of time during a procedure.
3. Resins should be refrigerated for long-term storage to maintain shelf life.
4. Light cured PQ1 is 15-25µm in thickness. Therefore, thin and air dry only before seating veneer with luting resin. For chemical cure or dual cure luting use the PermaFlo DC system.
C. General Precautions:
1. Replace luer lock caps when syringes are not in use.
2. When dispensing intraorally, never force syringe plungers. If there is excess resistance to flow, replace tip with a new one and verify flow before proceeding. Use only recommended tips with each material.
3. Tips should not be re-used so as to prevent cross-contamination.
4. To optimize bond strength, use oil-free, moisture-free air.

Quick Reference

STEPS:

BOND TO:	1 Surface Preparation	2 Apply Bonding Agent	3 Restore/Lute
A Dentin/ Enamel	**20-Second Etch** on enamel and dentin with Ultra-Etch®. Rinse firmly with air/ water spray. Blot or blow all visible water off but do not desiccate. A-1	Apply a puddled coat of PQ1 to entire preparation. Rub dentin with moderate pressure for 10 seconds, avoid rubbing enamel. Thin/dry with a focused air stream with ½ pressure 10mm for 10 seconds to remove alcohol solvent. Surface should look glossy, not puddled or milky. Light cure 20 seconds if using a standard halogen light, 10 seconds if using a high energy light, such as VALO. A-2	**Restore.** Light cure composite of choice. For amalgam, We recommend PQ Amalgam for highest strength amalgam bonding. A-3
B Porcelain Repair:	Freshen with diamond bur. Apply Ultradent® Porcelain Etch 90 seconds. Rinse and dry. Apply Ultradent® Silane 60 seconds, then air dry. B-1	SAME AS in A-2. B-2	SAME AS in A-3. If metal is exposed, address as in Db1. B-3
C Composite	Freshen with bur or microabrade. Rinse and dry thoroughly. Etch 10 seconds with Ultra-Etch to clean C-1	SAME AS in A-2. C-2	SAME AS in A-3. C-3

D ***Metal: a. Amalgam***	SAME AS in A-1. We recommend PQ Amalgam for highest strength in amalgam bonding.	SAME AS in A-2.	Pack amalgam onto cured bonding agent.
b. Cast	Best: Microabrade. Rinse and dry. Freshen with diamond bur if microabrasion is not available. Etch 10 seconds to clean. Rinse and air dry. D-1	SAME AS in A-2. D-2	SAME AS in A-3. First block out dark metals with a visually opaque, flowable composite such as PermaFlo Opaquer(s) and restore. D-3
Luting for Veneers	Same as in A-1	Same as in A-2.	Apply PermaShade LC to silanated veneer; seat to place, wipe off excess resin. Apply DeOx to prevent oxygen inhibition on margins and light cure.

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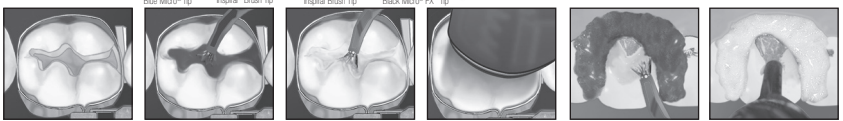
EN

Dentin/Enamel or Composite: Rubber dam placed, cavity prep completed.	Apply Ultra-Etch® 20 seconds. Air/water wash 5 seconds. Remove surface water, do not desiccate.	Rub PQ1 primer/bonding resin 10 seconds into dentin. Gently air thin 10 seconds. Surface should appear glossy and clear.	Light cure 20 seconds if using a standard curing light or 10 seconds if using a high energy light. Restore tooth using Amelogen® Plus, amalgam or other restorative material.	Porcelain: Apply Ultradent® Porcelain Etch 90 seconds. Rinse and dry. Apply Ultradent® Silane; air dry after 60 seconds. Bond (Figs. 3 & 4).	Cast metal: Microabrade (tin plate precious metals). Rinse and dry. Bond (Figs. 3 & 4).

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EN

General Precautions:

1. For Professional use only.
2. Review instructions, precautions, and MSDS before beginning treatment. Use only as directed.
3. Keep products out of heat/sunlight.
4. Avoid skin exposure to resins.
5. Isolate strong chemicals to area of treatment.
6. Confirm that patient has no known allergies to treatment materials.
7. Test flow of materials from syringe and tip before using intraorally.
8. Never force syringe plungers.
9. Use Ultradent® Syringe Covers and/or clean and disinfect syringes between patients.
10. Use standard curing lights with all light polymerized/activated materials.

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

PQ1

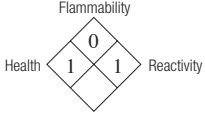For product MSDS please see our website: www.ultradent.com

Ultra-Etch

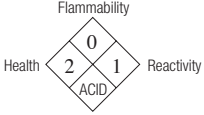For product MSDS please see our website: www.ultradent.com

HAZARD RATING

4 = Severe
3 = Serious
2 = Moderate
1 = Slight
0 = Minimal

Ultradent syringes have an expiration date stamped on the side of the syringe consisting of one letter and three numbers. The letter is a lot number used for manufacturing purposes and the three numbers are the expiration date. The first two numbers are the month, and the third number is the last number of the year.

U.S. federal law restricts the sale of this device by or on the order of a dentist.
Keep out of reach of children.

For immediate reorder and/or complete descriptions of Ultradent's product line, refer to Ultradent's catalog or call Toll Free **1-800-552-5512**. Outside U.S. call (801) 572-4200 or visit www.ultradent.com.

© Copyright 2011 Ultradent Products, Inc. All Rights Reserved.

ULTRADENT PRODUCTS, INC.

Manufactured by Ultradent Products, Inc.
505 West 10200 South, South Jordan, Utah 84095, USA Made in USA

63472-U.9 072613