EndoREZ®
Obturation Kit

EN

Description:
EndoREZ is a biocompatible, methacrylate-based root canal sealer/filler. It is a two-part, mixed dual cure material with optional EndoREZ Accelerator. EndoREZ has hydrophilic characteristics providing penetration into dentinal tubules. The radiopacity is equivalent to gutta-percha, providing simplified radiographic interpretation. EndoREZ does not compromise the function of dentin bonding agents or luting resin polymerization.

Indications for Use:
EndoREZ is designed to be used with EndoREZ Points and/or gutta percha for the filling of cleaned and shaped root canals. EndoREZ, in conjunction with a master cone and accessory cone, provides optimum sealing. Although Endo-REZ Points are recommended, EndoREZ may be used with all conventional endodontic obturation techniques.

Directions for Use:
A. ENDOREZ OBTURATION PROCEDURE WITHOUT OPTIONAL ENDOREZ ACCELERATOR
1. Instrument, clean and disinfect canals thoroughly.
CAUTION: DO NOT USE PEROXIDE, IODINE OR OTHER FREE RADICAL CONTAINING IRRIGANTS OR LUBRICANTS WITH ENDOREZ. IT WILL INHIBIT THE POLYMERIZATION OF ENDOREZ. CHECK PRODUCT INFORMATION BEFORE USING.
2. Evacuate liquid from canals for 2 seconds using Luer Vacuum Adapter with Capillary Tip.
3. Dry canal with paper points for 1-2 seconds each (point should be damp 1-3mm at tip). Canal should be damp, not desiccated, prior to obturating with EndoREZ.
CAUTION: DO NOT desiccate canal as this will cause less tubule penetration of sealer.
4. Fit EndoREZ Point master cone or traditional gutta percha point master cone to length, with tug back, for each canal.
5. Verify fit radiographically. Remove gutta percha from canal and set aside.
6. Remove EndoREZ cap from dual-barrel syringe by twisting counter clockwise. Discard.
7. Attach mixing tip to the dual-barrel syringe by lining up internal cap stems with syringe orifices. Twist clockwise to lock.
CAUTION: Dual-barrel syringe and mixing tip are NOT for intraoral use.
8. Remove small cap from mixing tip and express a small amount onto pad to verify flow.
9. Firmly attach appropriate length NaviTip to Skini syringe.
10. Insert mixing tip into back-end of Skini syringe and dispense enough EndoREZ into syringe to fill all canals.
11. Insert plunger and express air until material begins to flow from NaviTip.
12. Verify flow prior to applying intraorally. If resistance is met, replace tip and re-check.
13. Insert NaviTip into canal to within 2-4mm of apex or shorter if NaviTip binds in canal.
WARNING: To avoid expressing EndoREZ beyond the apex, DO NOT force NaviTip or material into tight canals.
14. Express EndoREZ with light pressure into canal while withdrawing tip. Keep NaviTip orifice buried in material while expressing EndoREZ and withdrawing tip.
15. Fill canal space to canal orifice ONLY.
16. Slowly insert the master cone (Step 5) and seat to working length. Harpoon in accessory cones to fill space as needed.
17. Light cure EndoREZ for 40 seconds.
NOTE: Initial surface polymerization with curing light (without Accelerator) is less than 0.3mm thick and aids in immediate restoration. Use caution not to disrupt surface seal when restoring. EndoREZ is fully polymerized in 20-30 minutes.
18. Trim excess gutta percha with a very hot instrument.
19. Restore.
CAUTION: Do not remove mixing tip until next use. Disinfect EndoREZ dual-barrel syringe and mixing tip with intermediate level disinfectant.

B. ENDOREZ OBTURATION PROCEDURE WITH OPTIONAL ENDOREZ ACCELERATOR
1. Follow steps A.1 – 16 above placing ONLY the master cone.
CAUTION: EndoREZ Accelerator will cause rapid cure of the EndoREZ resin and is intended for use only with accessory cones. Seat master cone and verify apical seal with radiograph PRIOR to inserting accessory cones with Accelerator.
2. Remove cap from EndoREZ Accelerator vial.
3. Dip a #25 or smaller accessory cone to its full length into the vial.
4. Insert accessory cone into filled canal as far as it will go. Place a minimum of 2-3 dipped accessory cones. Avoid excessive movement or pumping action while placing cones.
5. For larger canals, dip and insert additional accessory cones.
6. Re-cap Accelerator vial between canal uses. One vial is enough for approximately 3-5 canals.
7. Continue with steps A. 18 & 19 above.
8. Discard EndoREZ Accelerator vial after each patient.
NOTE: The maximum working time of EndoREZ with Accelerator is 3 minutes and will be fully set in approximately 5 minutes in canal.


C. ENDOREZ OBTURATION PROCEDURE WITH ACCELERATOR PLACED ON MASTER CONE
1. Follow steps A.1-15 above.
2. Dip pre-fitted master cone into accelerator vial slowly and passively and seat the master cone into the canal to working length.
3. Additional #25 accessory cones can also be dipped and placed passively into the canal.
4. Follow B.5-8 above.
NOTE: Once accelerator dipped master cone has been placed into the tooth polymerization will begin very quickly. This necessitates the clinician to be aware of root length and gutta percha position prior to initiation of polymerization.

Approved chemicals for use with EndoREZ:
Calcium Hydroxide, Citric Acid, EDTA, Chlorhexidine, Metacresylacetate, Formocresol, local anesthetic, saline, H20, Benzalkonium Chloride or Sodium Hypochlorite.

EndoREZ Warnings and Precautions:
1. Carefully read and understand all instructions before using. See individual product inserts for other precautions of products included in this kit.
2. Use only approved chemicals for irrigation, disinfectants, and lubricants. Other chemicals could inhibit set. DO NOT use iodoform gutta percha or irrigants and lubricants containing peroxide, iodine, or other free radicals when using EndoREZ. CHECK MANUFACTURER'S INFORMATION BEFORE USING. Thoroughly rinse ALL irrigants from canal with bacteria free water or local anesthetic as instructed above before placing EndoREZ in the canal.
3. Do not use on patients with known sensitivities to resins, including methacrylates.
4. The Dual Syringe and mixing tip are not intended for intraoral use.
5. Do not use plastic tips or Capillary tips to deliver EndoREZ into canal.
6. To avoid expelling EndoREZ past apex, keep delivery tip at least 2mm short of apical opening and loose in the canal. Do not "pump" delivery tip or gutta percha while expressing material into the canal.
7. To avoid clogging, use only recommended tips to deliver EndoREZ to the canal. Verify flow prior to applying intraorally. Do not overfill. While expressing EndoREZ and withdrawing tip, Clinician should note resin coming upward around the tip orifice.
8. Avoid cross-contamination of dual barrels. DO NOT attach and re-attach the same mixing tip.
9. Verify flow prior to applying intraorally. If resistance is met, replace tip and re-check.
10. Properly dispose of tip and empty syringe after use.
11. The EndoREZ mixing tip acts as the syringe cap between uses. Disinfect mixing tip and EndoREZ syringe with an intermediate level disinfectant after each patient. Change mixing tip BEFORE each use.
12. EndoREZ will not set as hard as standard resin composite.
EndoREZ can be removed for re-treatment or post placement.
EndoREZ will not set on mixing pad when exposed to air.

EndoREZ Accelerator Warnings and Precautions:
1. Carefully read and understand all instructions before using.
2. EndoREZ Accelerator vials are intended for single patient use only. Discard after each patient.
3. EndoREZ Accelerator is dispensed in spill-proof vials. However, use caution to keep vial upright while in use.
4. If gloves become contaminated with EndoREZ Accelerator, avoid contact with soft tissue and skin.
5. Avoid contact with soft tissue. If contact occurs, rinse with copious amounts of water.
6. Do not allow EndoREZ or Accelerator to be ingested.
7. Do not use EndoREZ Accelerator on a master come. Deliver to canal with accessory cones only.

NOTE: For MSDS and other information about using EndoREZ, EndoREZ Points (resin-coated gutta percha), EndoREZ Accelerator or other related products please go to www.ultradent.com.

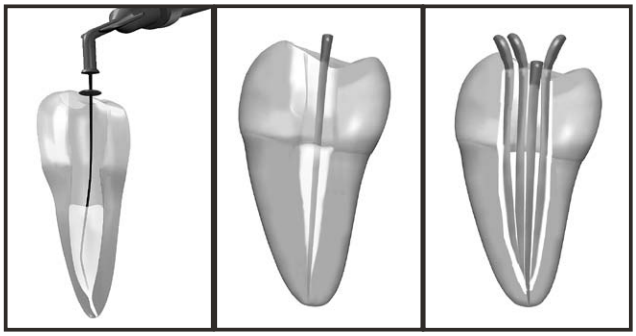EN Place NaviTip within 1-2mm of apex. Deliver short of apex with NaviTip loose in canal. Slowly insert master cone, then insert accessory cones as desired.

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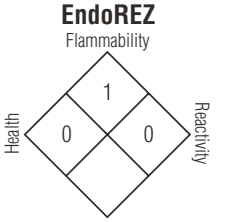
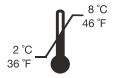
X026 = LOT
X026 = MM
X026 = YYYY

For immediate reorder and/or complete descriptions of Ultradent's product line, refer to Ultradent's catalog or call Toll Free1-800-552-5512. Outside U.S. call (801) 572-4200 or visit www.ultradent.com.

EndoREZ®
Obturation Kit

Manufactured by:
Ultradent Products Inc
505 West 10200 South
South Jordan, Utah 84095 USA
Made in USA

Refrigeration required
For professional use only.
Keep out of reach of children

For product SDS please see our website: www.ultradent.com

ULTRADENT PRODUCTS, INC.

© Copyright 2011 Ultradent Products, Inc. All Rights Reserved. 11595.7 012814