MTA FILLAPEX

PORTUGUÊS

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ENGLISH

MTA FILLAPEX is a root canal sealer based on MTA (mineral trioxide aggregate). Its composition provides the following qualities:

- Long term sealing capacity.
- High radiopacity.
- Promotes the deposition of hard tissue at the root apex and perforation sites.
- Low expansion during setting.
- Low solubility in contact with tissue fluids.
- Excellent viscosity for insertion.
- Does not stain the tooth.
- Eugenol-free material; does not interfere with the setting of resin cements.

COMPOSITION

Salicylate resin, diluting resin, natural resin, bismuth trioxide, nanoparticulated silica, MTA, pigments

INDICATION

For the final obturation of permanent teeth root canals, in combination with gutta-percha or silver points.

CONTRAINDICATION

In patients with hypersensitivity against the resins or other components of the product.

WARNINGS

- MTA FILLAPEX contains resins which may sensitize susceptible individuals. Do not use it in patients allergic to the resins or other components of the product.
- Avoid contact with eyes or skin. If it occurs, rinse immediately with water.
- Avoid contact with oral mucosa. If it occurs, rinse with water and prevent swallowing of product. In case any sensitivity persists, seek medical attention promptly.

PRECAUTIONS

- During the manipulation of this product, it is recommended the use of eye protection and appropriate clothes and gloves, for the patient and for the professionals.
- For better product conservation do not stock in refrigerator.

STORAGE


- Keep in dry and cool place, away from humidity and sources of heat.
- Keep the syringe tightly capped.

DIRECTIONS FOR USE

1. **Root canal preparation:** Prior to insertion of MTA-FILLAPEX, the root canal should be prepared, cleaned and dried, according to the selected endodontic technique.
2. **Mixing:** Performed by the self-mixing tip attached to the syringe. Use the cement immediately after mixing and dispensing by the syringe.

Warning: the self-mixing tip must be discarded after use.

3. **Insertion:** Use MTA FILLAPEX with gutta-percha or silver points, coating them with a thin layer of the cement and filling the root canal according to the selected technique. MTA FILLAPEX can be applied in the root canal with a Lentulo spiral or, directly, with an applicator tip adapted to the self-mixing tip.
4. **Removal of the root canal filling:** Use the conventional techniques for the removal of gutta-percha fillings.
5. **Working time:** 30 minutes.
6. **Setting time:** Minimum of 120 minutes (two hours).

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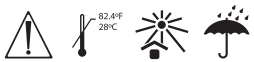ATTENTION: This product must be used according to the instructions described in this manual. The manufacturer is not responsible for failure or damage caused by incorrect handling or use.

ATENDIMENTO AO CONSUMIDOR
ATENCIÓN AL CONSUMIDOR
CUSTOMER SERVICE
+55 (43) 2101-3200
sac@angelus.ind.br
www.angelus.ind.br

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